




Deluxe CR Scanner Cart

Ergonomically designed to offer a practical workspace for the CR user as well as providing the mobility your site demands.

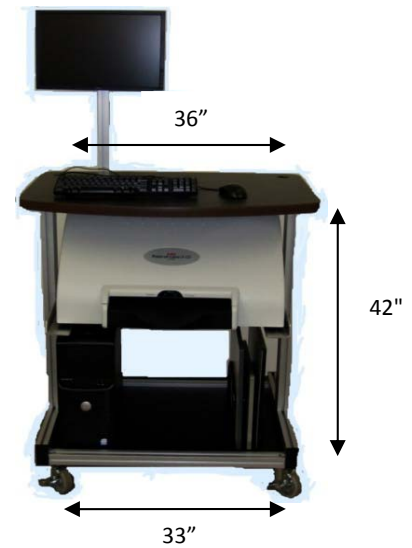

TI-BA's Deluxe CR Scanner Cart is designed specifically for Kodak Point-of-Care and Orex ACL CR Systems.

Features a sliding vertical/horizontal monitor rail system. LCD can be positioned at varying heights on either side of the cart. Extruded wire channels in the aluminum construction provide room to route cables.

Plastic snap-in covers secure wires in place to hide them, which adds to its aesthetics.

Table top measures 36" across and 24" deep to provide ample workspace. CR scanner rests on supports which allow access to the service panel underneath, without having to lift or move the unit. Lower shelf is tall enough to accommodate a server and UPS with built-in cassette racks to store and protect your valuable cassettes. Unit rolls on heavy-duty casters that can be locked once positioned.