

Improved!

Kodak Min-R S Mammography Film

Quality mammography imaging for challenged processing environments

KODAK MIN-R S mammography film is our most robust film designed for optimal performance with shallow-tank and mixed processing environments. This green-sensitive, high-contrast film features stable processing tolerance even in non-optimized environments. The result is mammography image quality you can trust.



Fits your environment

KODAK MIN-R S film requires no adjustment to equipment or procedures in most cases. Compatible with all KODAK MIN-R screens, this film:

- ▶ Provides a 100-speed system with MIN-R screens and a 150-speed system with MIN-R 2000 screens.
- ▶ Provides good, consistent results with shallow-tank and shared processing equipment.
- ▶ Can be used with orthochromatic mammographic screens from other manufacturers.
- ▶ Provides optimal results when used with the complete KODAK system including processing solutions, screens, cassettes, and processors.

Fine detail perception

Images are sharp and provide viewing ease of diagnostic information. Perception of microcalcifications and other minute details is enhanced by the higher contrast and sharpness of the film.

New innovative design


Now with a pleasing blue tone, KODAK MIN-R S film is an innovative dual-emulsion film designed for use with a single intensifying screen. A proprietary technology employs two thin emulsions working in tandem to yield a superb image—while providing greater processing tolerance than a thicker single emulsion.

- ▶ The primary emulsion images breast parenchyma, providing detail and now even higher contrast.
- ▶ The secondary emulsion images the breast periphery and upper density.


A name you can trust

Carestream Health is a leader in mammography imaging because of our clinical expertise and commitment to your satisfaction. We combine state-of-the-art film technology with the right processing chemicals and equipment to provide the medical images you need for high-quality patient care.

Dual-emulsion technology

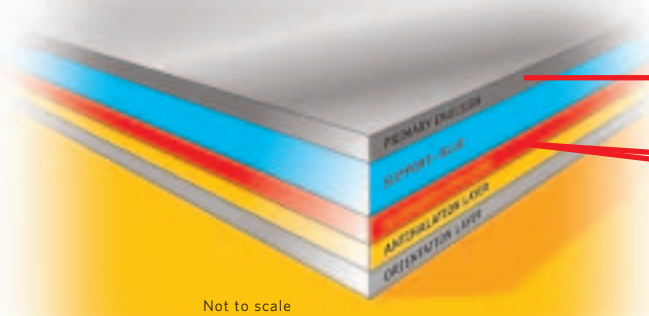

Not to scale

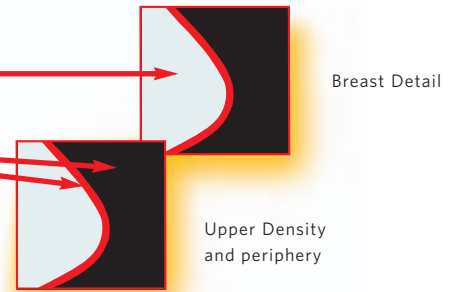

Ordering information

PRODUCTS	CAT NUMBER
KODAK MIN-R S film 24 x 30 cm	173 2726
KODAK MIN-R S film 18 x 24 cm	151 3381

KODAK MIN-R 2 cassette with:

KODAK MIN-R 2000 screen 24 x 30 cm	899 9492
KODAK MIN-R 2000 screen 18 x 24 cm	810 4101
KODAK MIN-R 2190 screen 24 x 30 cm	891 0853
KODAK MIN-R 2190 screen 18 x 24 cm	126 0579
KODAK MIN-R screen 24 x 30 cm	709 1614
KODAK MIN-R screen 18 x 24 cm	708 7984

KODAK MIN-R cassette with:

KODAK MIN-R 2000 screen 24 x 30 cm	193 0247
KODAK MIN-R 2000 screen 18 x 24 cm	878 4449
KODAK MIN-R 2190 screen 24 x 30 cm	874 3171
KODAK MIN-R 2190 screen 18 x 24 cm	865 0426
KODAK MIN-R screen 24 x 30 cm	128 4413
KODAK MIN-R screen 18 x 24 cm	128 4397

Screen-film characteristics

KODAK Screen	KODAK Film	Relative Processing Cycle	Relative Speed	Contrast		D-Max
				RP	EXII	
MIN-R 2000	MIN-R S	Standard	150	3.80	4.0	≥ 4.0
MIN-R 2190	MIN-R S	Standard	190	3.80	4.0	≥ 4.0
MIN-R	MIN-R S	Standard	100	3.80	4.0	≥ 4.0

Part of an unbeatable system

Use the MIN-R S film system with Kodak processing chemicals and equipment for mammography to help your facility maintain the highest standards of performance and care.

Kodak

Licensed Product

Carestream Health, Inc.
150 Verona Street
Rochester, NY 14608

More information

For details on KODAK MIN-R S mammography film, or to order, contact your Carestream Health representative or distributor, or call 1-877-865-6325, ext. 626.

→ www.carestreamhealth.com

Carestream 
HEALTH

The Innovation Powering **Kodak** Health Products