

Detector Storage Bin

Detector Locker
Anti-Theft Lock Function

Battery Box

Built-in Charger
Universal Connector
Easy Magnetic Connector
Contamination Free



Analog Type for CR or Film Users



Control Console
10.1" Touch Screen

Options

- Fixed Column
- DAP(Dose Area Product) Meter
- AI Assisted Diagnostic Imaging
- Scatter Reduction
- Bone Suppression

Specification

High Frequency X-ray Generator

Model	PROMO
Input Power	100 - 240 V~ (Free Voltage Input), 3kVA, 50/60 Hz
Rated Output	40 kW
mA Range	Max. 500 mA
kV Range	40 - 125 kV

promo

Mobile X-ray System

*Experience Positioning Flexibility,
Superior Performance*



DRGEM
Your Best Healthcare

© DRGEM Corporation - All right reserved.
Reproduction in whole or in part is prohibited without the written consent of the copyright-holder.
Specifications and features are subject to change without notice.
All changes will be in compliance with regulations governing manufacture of medical equipment.

Head Office/R&D 7F, E-B/D, Gwangmyeong Techno-Park, 60 Haan-ro, Gwangmyeong-si, Gyeonggi-do, Korea (14322)
Gimcheon Factory 213, Saneopdanji-Ro, Eomo-Myeon, Gimcheon-Si, Gyeongsangbuk-do, Korea (39536)
Gumi Factory 116-59, Sanho-daero, Gumi-si, Gyeongsangbuk-do, Korea (39377)

DRGEM
Service Community

24 hours 7 days
Customer Support
through Salesforce
Service Cloud CRM

DRGEM
Your Best Healthcare

Positioning Flexibility and Seamless Workflow

PROMO

The PROMO

PROMO enables positioning flexibility, easy maneuverability, and seamless workflow, ensuring high-quality images every time

High-Quality Imaging at Point of Care

PROMO delivers exceptional imaging services directly to the bedside, including in the ER, OR, ICU, and patient rooms. Its flexible 180-degree rotating column and maximum X-ray focal spot height facilitate easy positioning, quicker exams, and a faster workflow

Advanced Imaging Technology

Equipped with a powerful 40kW capacitor generator, dual-focal spot tube, and scatter reduction image processing features, PROMO optimizes image quality, ensuring clear and precise diagnostic images. Its on-device AI assists in the immediate detection of suspected chest abnormalities indicative of major pulmonary diseases from X-ray images


User-Friendly Interface

The push-down-and-go driving system, on-site image verification on a large 21.5-inch battery-powered touch monitor, remote control, state indicators, multi-purpose storage, built-in detector charger, and intuitive user interface allow for fast examinations and improved operability, significantly reducing training time

Digital Efficiency and Connectivity

PROMO offers the speed, efficiency, and consistency of digital imaging. Operators can easily access the Hospital Information System (HIS) from any location, enabling seamless communication with the Radiology Information System (RIS) and Picture Archiving and Communication System (PACS)

Rotating Column

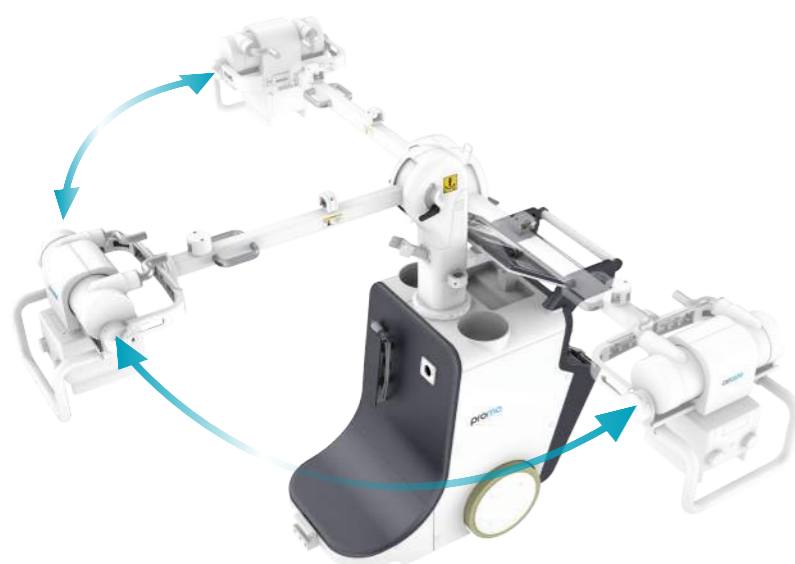

Highest Focal Spot



Rotating Column & Articulated Arm
Maximum Flexibility and Reach

Remote Control & Receiver

Collimator Lamp, Ready, and Exposure control



Touch Monitor

21.5" Touch Monitor
Built-in Computer
PC Alive Feature

Handle Bar & Brake

Push Down and Go
for Easy Driving


Exposure Switch

Detector Tray

Built-in Charger
Mechanical Locker
Two Detector Storage

Kickstand

Quick Lift Over Doorways
and Elevator Thresholds

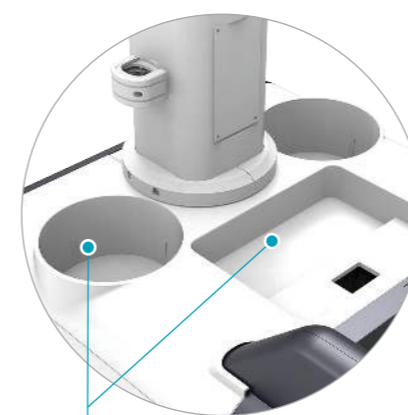

Multi Purpose Storage

Isopropyl Alcohol Wipes
Small Size Sterile Gloves
Lead Markers, Tapes and etc.

Additional Interface

USB / LAN / Exposure Switch / Connector



State Indicators

LED State Indicators
/ Power Switch

